

Your best solution for the removal of thrombus

The Escape Aspiration Catheter is a single-use, rapid exchange, dual lumen catheter for the removal of thrombus in the central and peripheral circulation system.

■ *Optimised tip design*

Smooth tip design with integrated radiopaque tip marker to increase deliverability

■ *Excellent pushability and kink resistance*

Innovative shaft and stylet support to avoid kinking

■ *Optimal trackability*

Contour™ hydrophilic coating on 25cm long outer lumen

■ *Large extraction lumen*


Uniformly large extraction area for increased aspiration speed

ESCAPE


Aspiration Catheter

Technical Data	Escape (6F)	Market Leader
Aspiration Speed	1.6cc/s	1.1cc/s
Extraction lumen area	0.95mm ²	0.81mm ²
Outer diameter (O.D)	0.067" (1.70mm)	0.066" (1.68mm)
Proximal shaft push/kink resistance	Reinforced PEEK and stylet	Braided shaft

- 1** Aspiration Syringe
- 2** Filter
- 3** Needle
- 4** Stop cock

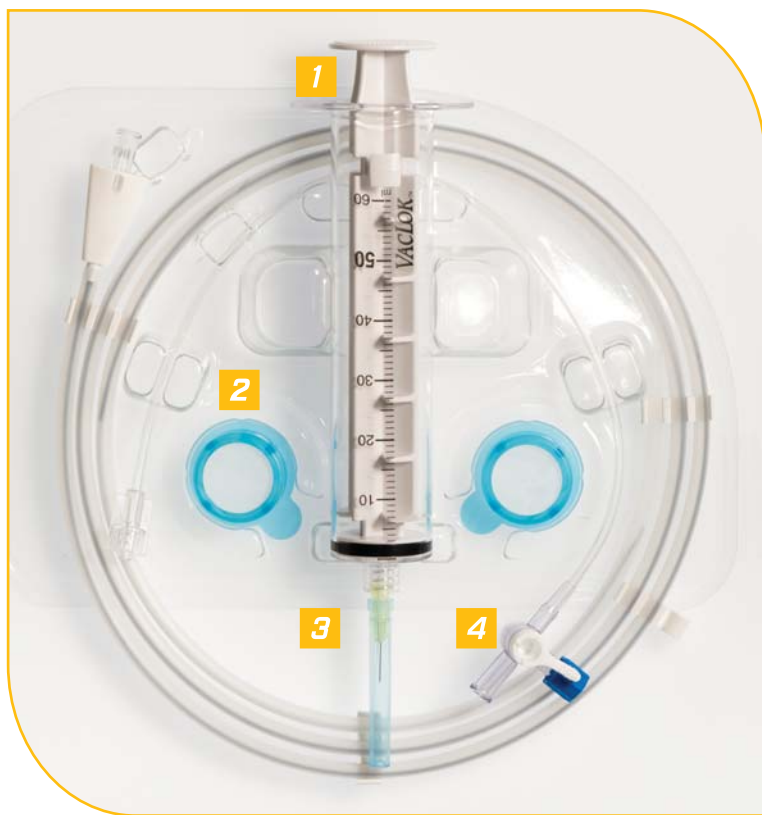

- Catheter length: 145cm
- Rapid exchange length: 25cm

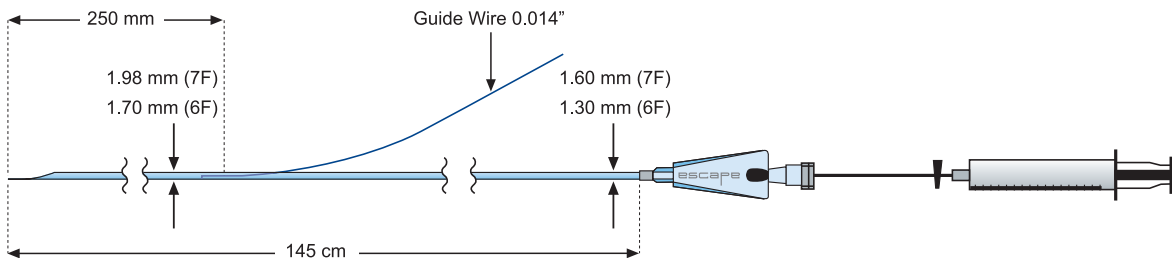

Please see the instructions for use for a complete listing of the indications, contra-indications, warnings and precautions.
 © 2011 Escape is a trademark of DISA Vascular
 Escape Aspiration Catheter Brochure Issue 02/01/12

www.disavascular.com



Bringing Quality to Life